

# Sonic

WHITENING TOOTHBRUSH

USER MANUAL


**b** beconfident®

# Congratulations on your quality purchase from Beconfident® – We hope you will love our Sonic toothbrush and toothpastes!


Before using this product, read this user manual carefully to learn about the features and the recommendations (the do's and don'ts) and how you can make the best use of it. **Note:** Charge the Beconfident sonic toothbrush completely prior to using it for the first time.

## Quick guide how to use the brush:

1. Press power button to turn on the brush
2. Press again in less than 5 seconds to toggle through brushing modes.
3. Long press power button for 2 seconds to turn off the brush after your favourite mode is selected.
4. Wet the bristles and apply a small amount of toothpaste.
5. Place the toothbrush bristles against the teeth and (do not press hard).

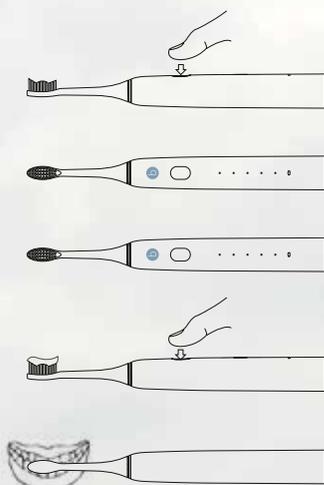

Register your purchase at [www.beconfident.info/sonicregistration](http://www.beconfident.info/sonicregistration) with purchase date and serial number for extended warranty and receiving Beconfident Sonic Newsletter for tips and special offers.

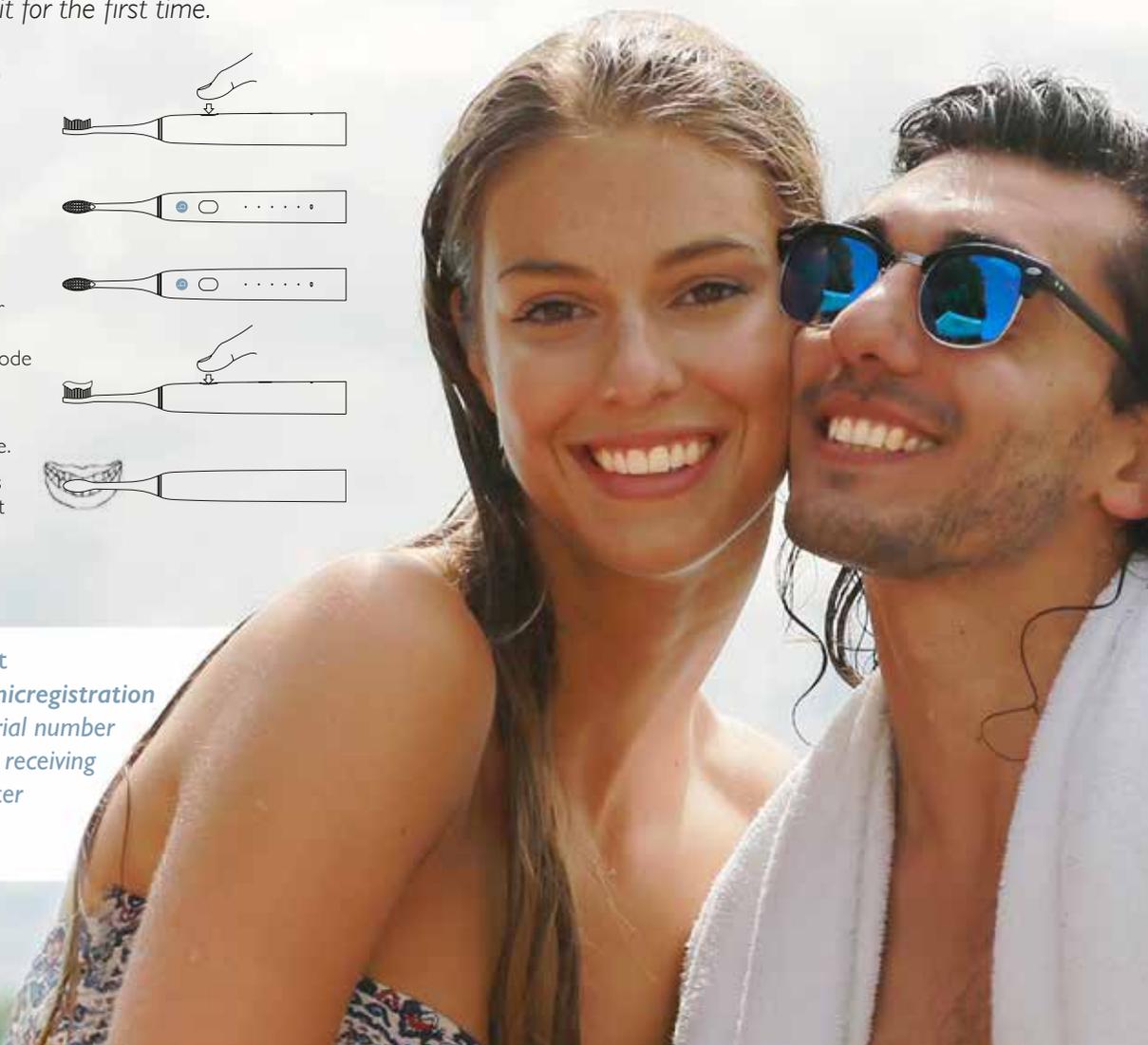

## Important safeguards

**Warning:** Read all instructions before use

1. Check the product carefully before use every time. Discontinue use if the product appears damaged in any way (brush head, handle, charging pad).
2. The charging pad cord cannot be replaced. If the cord is damaged, discard the charging pad and purchase a new one.
3. This product contains no user-serviceable parts. Do not disassemble, repair or modify the product on your own. Only certified professional technicians are qualified to fix any problems. Refer to "WARRANTY AND SUPPORT" if the product no longer works properly or needs repair.
4. Always place and store the product (including charging pad) far away from children or the mentally disabled.
5. Always store the product in a cool and dry place away from direct sunlight when you don't use it for an extended period of time.
6. To avoid damage to the product, do not place the brush head (including protective cap), handle, charging pad in the dishwasher for cleaning.
7. Replace Beconfident brush heads every 3 months to achieve optimal results. Use only Beconfident replacement brush heads.
8. Do not use the brush while charging.
9. The product is supplied with a rechargeable Li-ion battery. This battery is self-contained and not consumer replaceable. Battery must be disposed of properly when it no longer holds a charge.
10. Always place and store the USB plug of the charging pad at a safe distance from sources of water (i.e. sink or tub) although the charging pad main part is IPX-7 waterproof.
11. Do not put the toothbrush on the charging pad that has fallen into water. Unplug immediately.
12. Never use a charging pad with a damaged cord.
13. Do not use while warm/hot bathing or shower. Although the product is IPX7 waterproof, it cannot be effectively protected from hot water or steam.
14. This product is designed to clean your teeth, gums and tongue only. Use this product only for its intended use as described in this user manual or as recommended by your dental professional.
15. This product is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the electronic device by a person responsible for their safety.
16. Children under the age of 7 are not recommended to use this product.
17. This product is a personal care device and is not intended for use on multiple patients in a dental office or institution.
18. Use UL listed computer, USB hub or USB duplex receptacle when charging the brush handle.
19. In order to effectively charge your brush do not place your charging pad on a magnetic metal surface.
20. For connection to a supply not in your country, use an adapter of the proper configuration for the power outlet.
21. Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
22. Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after one week of use.
23. Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
24. This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.
25. Consult your physician if you have other medical concerns.
26. This brush is safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers).
27. It is normal for your gum to feel itchy, numb or even slightly bleeding after your first-time use. Your gum will get used to the sonic vibration in about two days and the discomfort will decrease in one week. Consult your physician if discomfort continues.

**Product:** Beconfident Sonic Whitening Toothbrush  
**Rated Voltage:** 3.7V  
**Rated Power:** 2.0W  
**Responsible Company:** Group Beconfident Europe AB  
**E-mail:** [info@beconfident.info](mailto:info@beconfident.info)  
**Product Standard:** GB4706.1-2005, GB 4706.59-2008

## Items in the box

- 2 Brush heads with hygienic transparent cap.\*
- Handle.
- Charging station (without adapter).
- User manual.

\* Please do not dispose of the transparent cap as it can be used to protect the bristles of the head in case of travels.


**Notes:** There is no adaptor included in the box. If needed, any smartphone or other electronics adaptor with output voltage 5V can work safely with the charging pad for charging the brush.




## Specifications

Model No	White: 650098, Black: 650598
Size	251x29x28.5mm
Frequency	31,000~40,000VPM
Brushing modes	Clean, Sensitive, Gum Care, White, Polish
Working voltage	3.3~4.2V
Battery.	2,000mAh Li-ion, last 90 days on single charge
First time charging	Charging takes about 30 hours to full charge
Low battery warning	Flashing red + vibrating motor at the moment of switching on
Power button	To toggle through modes press the power button twice in less than 5 seconds. To turn off your brush press the power button again after 5 seconds.
Indicator on charging	When charging: Solid red. When fully charged: Solid green for 30 sec and then off
Waterproof	IPX-7
Smart memory	The brushing mode of your last use will be remembered and activated in next use.
Brush head	Food-grade Dupont bristles with protective cap
Noise	<60dB
Weight	139g (with brush head)

## Product diagram

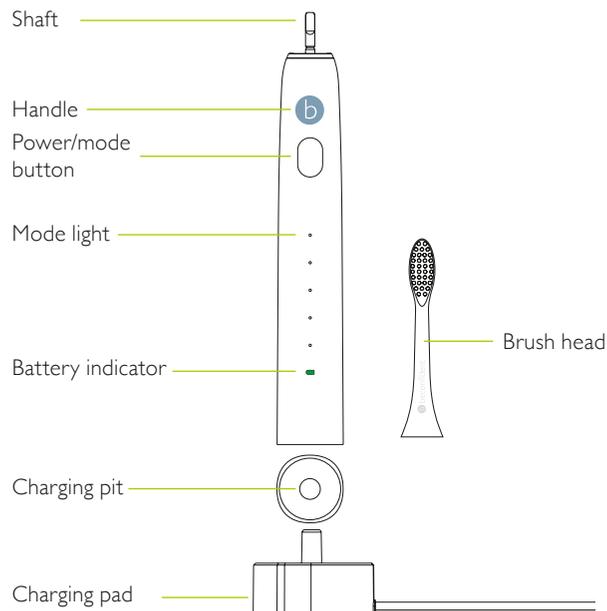


## Brushing modes

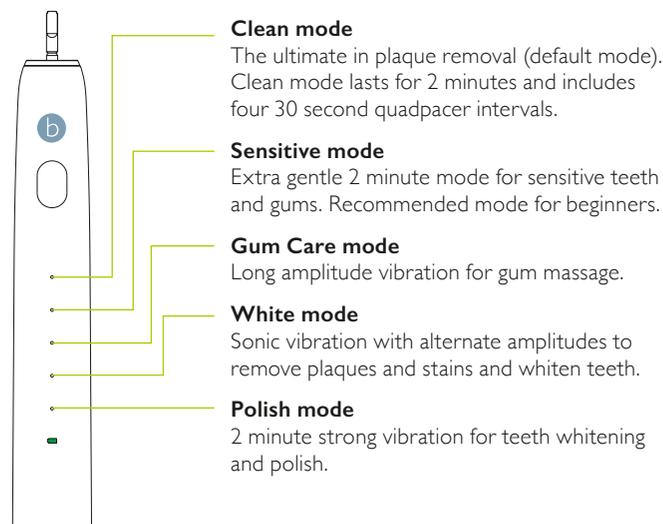

## How to install the brush


Align the brush head so the bristles point in the same direction as the front of the handle. Firmly press the brush head down onto the metal shaft until it stops.

It is normal to see a slight gap between the brush head and the handle.

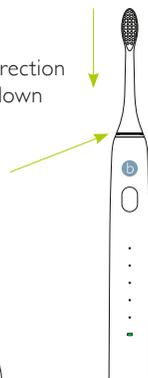


## Weekly maintenance

Rinse the top of the brush handle and the bottom of the toothbrush head once a week.

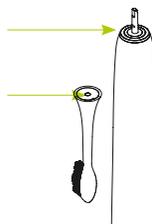

## How to use the brush

1. Press power button to turn on the brush.
2. Press again in less than 5 seconds to toggle through brushing modes.
3. Long press power button for 2 seconds to turn off the brush after your favourite mode is selected.
4. Wet the bristles and apply a small amount of toothpaste.
5. Place the toothbrush bristles against the teeth and (do not press hard).

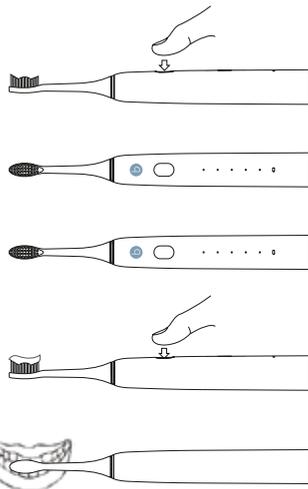



## How to turn off the brush

1. After you have completed the brushing, it will be off automatically.
2. Long press the power button for 2 seconds to turn off the brush.
3. Press the power button to turn off the brush after working in more than 5 seconds.

## How to brush your teeth

1. Every 30 second quadpacer interval reminds you to brush each section of your mouth.

**Notes:** To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the quadpacer feature.

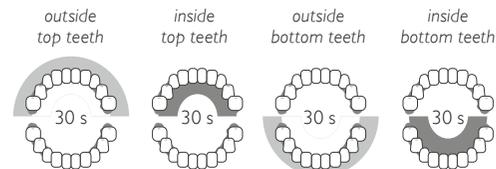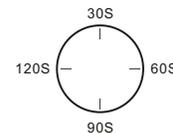

2. For optimal cleaning always place your brush at a slight angle (45 degrees) towards the gum.



3. Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion.




4. Pull out the brush head in the arrow direction for head replacement.



7. Rinse your brush head and handle with water after each use. Wipe the entire surface of the handle with a soft tissue or cloth.



5. Make sure you turn off the brush before taking it out of your mouth to prevent splashing saliva and toothpaste.

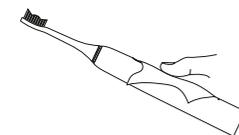


6. If you are not going to use your brush for an extended period of time, clean it and store it in a cool and dry place away from fire, heating objects or direct sunlight.



8. Charge the brush only when needed. Normally every 60-90 days.

*Power Indicator*  
When charging:  
Light is red

When fully charged:  
Light is green  
for 30 sec then off

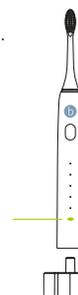

## Warranty and support

Beconfident warrants its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Beconfident's (including Beconfident retailers and dealers) expense provided that convincing proof of purchase in the qualifying period is provided. Beconfident recommends use of genuine Beconfident brush heads for optimal performance. Use of other brush heads may damage your brush and limit your warranty protection. Contact your local retailer and dealer or send an Email to: [info@beconfident.info](mailto:info@beconfident.info) for warranty service.

## Warranty exclusions

What is not covered under warranty:

- ✓ Brush heads.
- ✓ Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- ✓ Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- ✓ Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.
- ✓ The proof of purchase has been altered, defaced or removed.

All implied warranties, including warranties of merchantability and fitness for a particular purpose are limited in duration to the length of this warranty. Beconfident or any of its affiliated or subsidiary companies shall not be liable for any special, incidental, consequential or other damages based upon breach of warranty, breach of contract, negligence, tort, or any other legal theory. Such damages include, without limitation, loss of savings or revenue, loss of profit; loss of use; the claims of third parties including, without limitation, dentists and dental hygienists; and cost of any substitute equipment or services.

If the Beconfident product is required to be sent back to the retailer or dealer for repair or replacement, the following must be furnished to have the product serviced and returned:

1. For warranty return include some form of proof-of-purchase, such as a mechanical reproduction or carbon of a sales receipt. Make sure the date of purchase and product model number are clearly readable. If the originals are sent, they cannot be returned.
2. Send the entire product after contacting customer service.
3. Pack product securely to prevent damage in transit. If possible, use the original packing material.

4. Enclose a description of what was happening with the product. Include a typed or clearly printed name and address of where the product is to be returned, with phone number (required for shipment). Please fill out below customer info form:

Date of Purchase: \_\_\_\_\_

Email: \_\_\_\_\_

Address: \_\_\_\_\_

Description of the Issue: \_\_\_\_\_

Name: \_\_\_\_\_

Telephone: \_\_\_\_\_

- ✓ This product is warranted for a period of 24 months from date of purchase.
- ✓ Brush Sonic – Beconfident!

## Daily whitens and cleans teeth carefully

- ✓ Sonic Whitening Toothbrush.
- ✓ Sonic Whitening toothbrush with 38000 oscillations per minute.
- ✓ 5 brushing modes - Clean - Sensitive - Gum Care - Whitening and Polish enamel.
- ✓ 2-different types brush heads included.
- ✓ Travel battery built in for 90 days when used 4 minutes daily.
- ✓ One button for power on/off and mode selection.
- ✓ 2 minute timer with 30 second interval pause.
- ✓ Prevents cavities, tartar & stains.
- ✓ Wireless charging mode via USB.
- ✓ 24 Month standard warranty.
- ✓ Sonic is IPX7 waterproof. CE, Rohs.FCC and FDA approved. Designed in USA and distributed by Beconfident worldwide.\* Extended ww warranty.

Beconfident® is a registered trademark of Group Beconfident Europe AB  
[www.beconfident.info](http://www.beconfident.info)



# BECONFIDENT MULTIFUNCTIONAL TOOTHPASTE



- Whitens teeth within two weeks
- Strengthens enamel
- Freshens breath
- Protects against plaque & cavities

**b** beconfident®

Group Beconfident Europe AB  
Tennisvägen 31  
SE-269 33 Båstad, Sweden

www.beconfident.info  
.com .co.uk .de .fr .nl .be .es .it .se .no .dk .fi